Model #TV865 | Modelo #TV865

TV LISTENER[®] K4

Rechargeable Wireless Headset
Auriculares Inalámbricos Recargables

Instruction Manual

Please save for future reference.

Manual de Instrucciones

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1-877-736-6760 / www.unisar.com

Warning

Please read all safety and operating instructions carefully and completely before using this product.

1. For your well-being, medical experts advise against playing any audio device at a high volume for continuous and extended play. If you hear a ringing in your ears, reduce the volume immediately or discontinue use.
2. To reduce the risk of fire or electrical shock, do not expose this product to water or moisture.
3. To prevent overheating, keep this product away from heat sources such as radiators, stoves, fireplaces, or direct sunlight.
4. Use only the AC adapter provided with the transmitter. Use of any other adapter may damage the unit and void the warranty.
5. Plug the AC adapter into standard household current only (120 volt AC outlet).
6. Avoid damage to the adapter cords. Route cords so that they are not likely to be walked on or pinched by items placed upon or against them.
7. Use only AAA NiMH rechargeable batteries with the headset.
8. Do **NOT** use alkaline batteries or you will severely damage the unit.
9. Dispose of batteries properly and safely.
10. Unplug the AC adapter and remove batteries during long periods of non-use.
11. Do **NOT** disassemble.

Thank you for purchasing the Unisar® TV Listener® K4. This infrared wireless headset system will enable you to listen to your television, stereo, computer, or any other audio device at a volume that is comfortable for you without disturbing others in the room. It is intended for use if you have difficulty hearing, if you simply want to watch television at a different volume than others in the room, or if you want to watch television while someone else is sleeping.

Please review the installation choices that are in this manual and determine which features your television or audio device has before setting up your TV Listener K4. Please save this manual for future reference.

If you have any questions, please contact our Customer Service Department at **1-877-736-6760** Mon-Fri from 9:00AM-5:00PM Eastern Time, or e-mail us at customerservice@unisar.com.

Table of Contents

How It Works / page 4

Setting Up the System / page 4

Operating the System / page 9

Tips for Using Your TV Listener K4 / page 11

Care and Cleaning / page 11


Troubleshooting / page 12

Technical Specifications / page 14

Limited Product Warranty / page 15

Also Available / page 16

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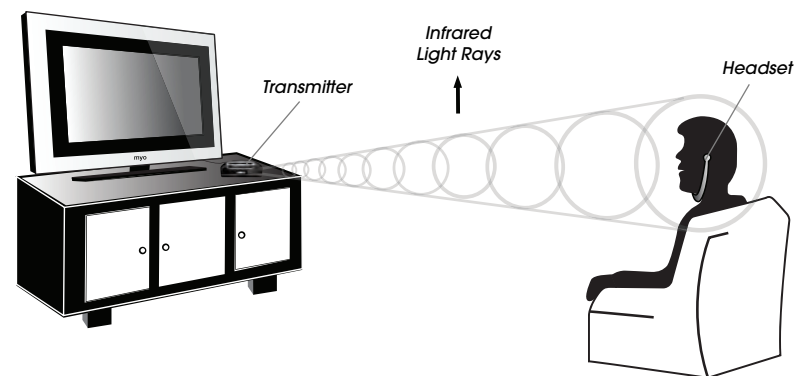How It Works

The TV Listener K4 works just like your remote control, sending an invisible infrared signal from your television or any other audio device to the headset. Unlike radio frequency devices, infrared is not subject to interference from cordless phones, wireless networks, microwaves, etc. You adjust the volume on the wireless, lightweight headset to suit your personal needs. Multiple headsets can be used with one transmitter.


Setting Up the System

Charge the batteries in the headset for 8 hours prior to first use.



TV Listener K4 Components

Check that you have all the following components before installation.

A Transmitter (connects to your TV or other audio device)

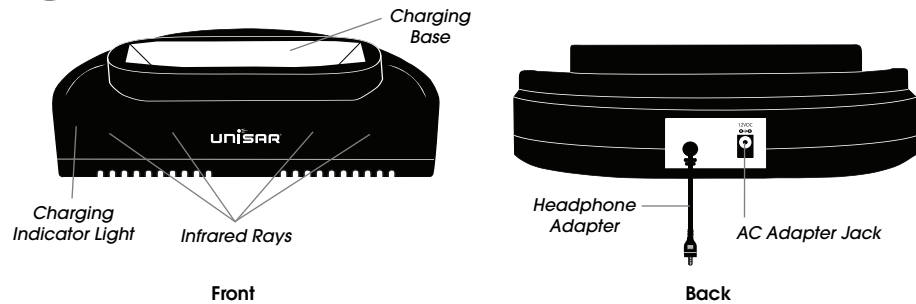B Wireless Headset

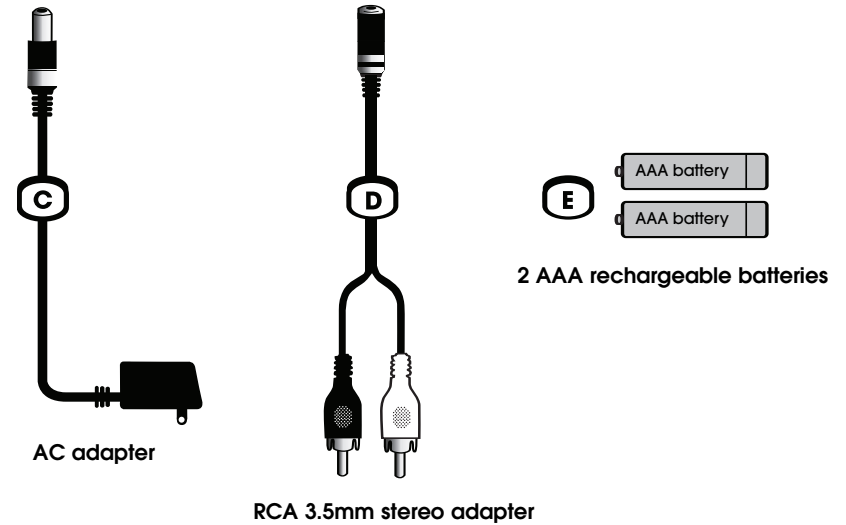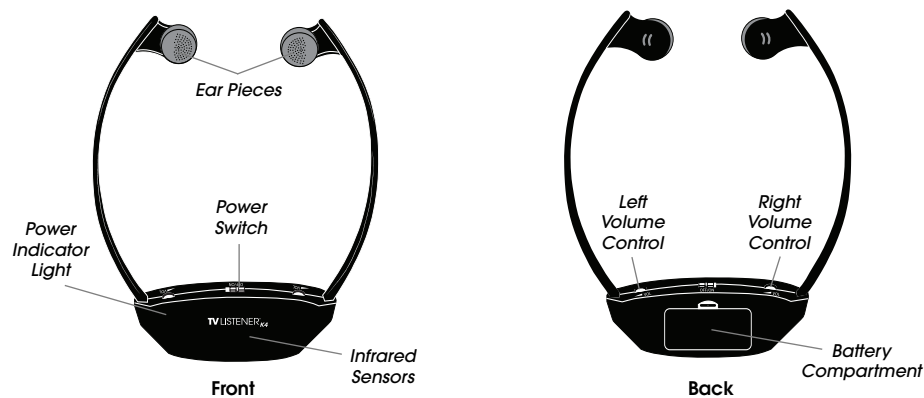
Installation of the Transmitter

The following are the two different ways to connect the TV Listener K4. Your installation will depend on the type of equipment you have.

1. **STEREO INSTALLATION** – for equipment with two audio output jacks (stereo televisions, DVD players, cable boxes). Refer to Pg. 6.
2. **SMALL HEADPHONE JACK INSTALLATION** – for equipment with a small headphone jack (most CD players, portable stereos, MP3 players). Refer to Pg. 7.

STEREO INSTALLATION

Fig. 1

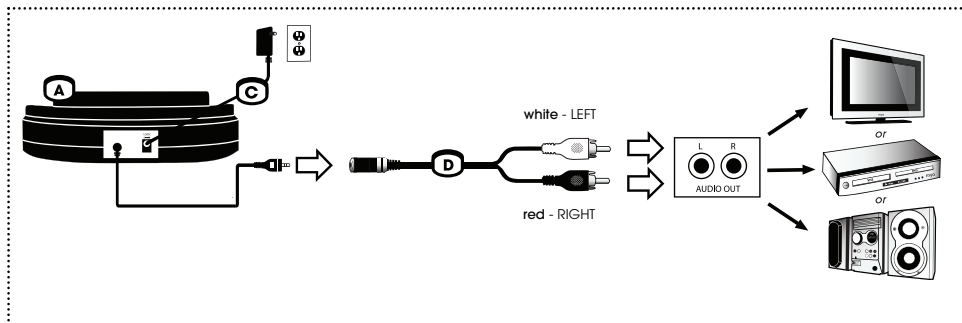Follow these instructions if your equipment has two audio output jacks (stereo televisions, DVD players, cable boxes).

1. Connect the audio cord on the back of the transmitter (**A**) to the RCA 3.5mm stereo adapter (**D**).
2. Plug the red end of the adapter into the right audio output jack and the white end into the left audio output jack on your equipment.


If you are connecting to your television, use the audio jacks located on the back. The jacks are labeled "AUDIO OUTPUT" (not just "AUDIO"). Do **NOT** use the audio jacks on the front as these are audio **input** jacks. Be sure to plug the adapter into the audio output and not the video output.

If you are connecting to a device that is connected to your television (DVD player, VCR, cable box, home theater system), this device must stay on when using the TV Listener. The VCR must also be set to the same channel as the television.

3. Plug the AC adapter (**C**) into the jack marked "12VDC" on the back of the transmitter and plug the adapter into an electrical outlet. For long periods of non-use, we suggest unplugging the AC adapter.

SMALL HEADPHONE JACK INSTALLATION

Fig. 2

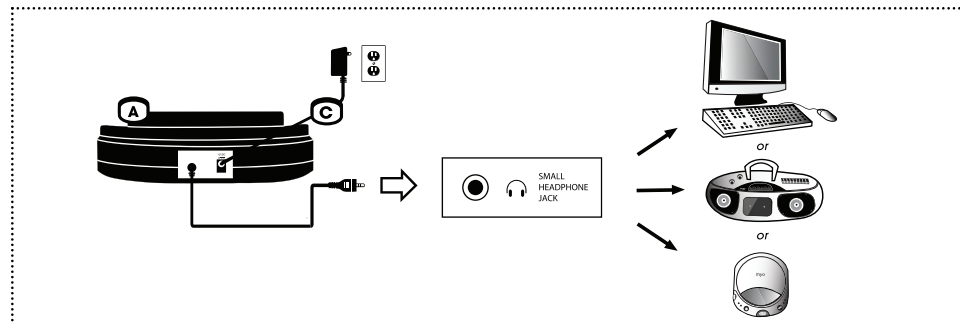Follow these instructions if your equipment has a small headphone jack (CD players, portable stereos, computers, MP3 players).

1. Connect the audio cord on the back of the transmitter (**A**) directly into the headphone jack on your equipment.
2. Plug the AC adapter (**C**) into the jack marked "12VDC" on the back of the transmitter and plug the adapter into an electrical outlet. For long periods of non-use, we suggest unplugging the AC adapter.

NOTE: If you connect the transmitter to your television using the headphone jack, the sound from your television will only be heard through your headset and will not be heard from your television speakers.

Headset Battery Installation

IMPORTANT: Charge the batteries for 8 hours prior to first use. You must only use rechargeable batteries (included). Do NOT use alkaline batteries or you will severely damage the unit.

The wireless headset (B) requires two AAA rechargeable batteries (included). To access the battery compartment, remove the battery cover on the back side of the headset (refer to Fig. 3). Insert two AAA rechargeable batteries matching positive (+) to positive (+) and negative (-) to negative (-) and replace the battery cover. Do NOT use alkaline batteries.

Charge the batteries for 8 hours prior to first use. Be sure the headset is turned off and the transmitter is plugged in to an electrical outlet. The transmitter will automatically turn on when the AC adapter is plugged in. Place the headset in the charging base of the transmitter (refer to Fig. 4). The red charging indicator light on the left side of the transmitter will illuminate. On a full charge, the battery life is approximately 10 hours. Recharge when the volume begins to weaken and the power indicator light on the headset begins to dim. It will take 6 hours to fully charge the headset after the initial charge. The charging indicator light will turn off once the headset is fully charged. If the charging indicator light does not turn on when you put the headset in the charging base, the headset is at least 70% charged. It will not charge until the headset is less than 70% charged. This helps to prolong the life of the batteries.

For long periods of non-use we suggest removing the batteries.

The batteries need to be replaced when they last much less than the expected time (approximately 10 hours) after proper charging. Be sure to replace **ONLY** with AAA NiMH rechargeable batteries. Do **NOT** use alkaline batteries.

Fig. 3

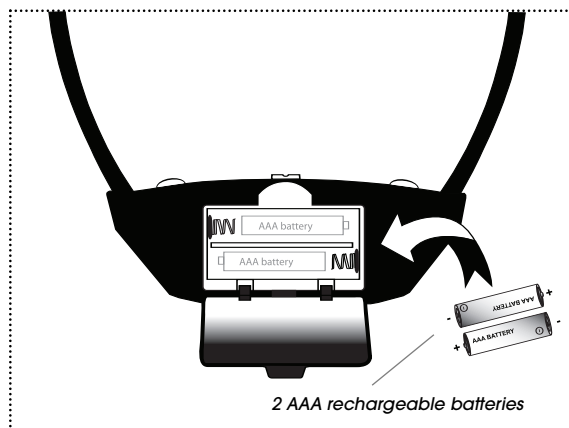
Fig. 4


Operating the System

1. Place the transmitter at or close to your eye level when you are sitting. There should be nothing obstructing the transmission of the infrared light rays from the transmitter to the receiver in your headset (refer to Fig. 5). The transmitter must be within 33 feet (10 meters) of the headset for the system to work properly.


NOTE: Avoid exposure to direct sunlight. Direct sunlight may interfere with the reception. It is also subject to interference from treadmills.

Fig. 5

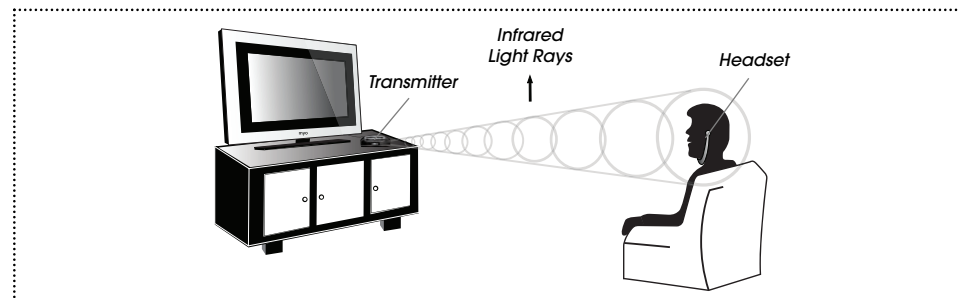2. The transmitter will automatically turn on when the AC adapter is plugged in. Four red lights on the front of the transmitter will illuminate.
3. Turn the headset on by sliding the power switch to the "ON" position. The red power indicator light will illuminate. To prevent possible damage to your ears, be sure to lower the volume on the headset before putting it on.
4. Put the headset on making sure the right ear piece (R) goes in your right ear and the left ear piece (L) goes in your left ear. The ear pieces do not go inside the ear. They rest just outside. The headset should hang under your chin with the "TV Listener K4" logo facing away from you (refer to Fig. 6).

Fig. 6

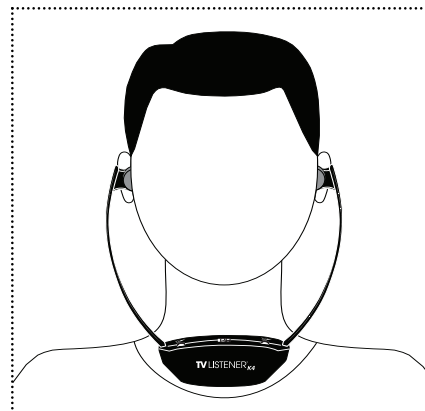5. Adjust the volume for each ear to a comfortable listening level using the volume control dial. If the headset is out of the range of the transmitter or the infrared sensors in the headset are being blocked, no sound will come through the headset.
6. Turn off the headset after each use and place it in the charging base to recharge. The transmitter will automatically stop charging once the headset is fully charged so the headset can remain in the base. We recommend unplugging the transmitter from the electrical outlet and removing the batteries from the headset during long periods of non-use.
7. The four red lights on the transmitter will turn off and go into standby mode if no audio signal is detected for three minutes (either you unplugged the transmitter from the TV or shut off the TV). Once the audio signal is resumed, the lights will turn back on. This feature helps to reduce power consumption.

Regardless of what type of equipment you have and which installation instructions you follow, when you connect the transmitter the sound will be heard either (1) through your headset and television speakers, or (2) just through your headset. The following instructions will apply to both types of installation.

1. **The sound from your television is heard from your headset and your television speakers.** If you want the sound to be heard only from the headset, you will need to mute the television sound. If you simply turn the volume down all the way, you may reduce the signal going to the television audio output jacks which will affect the sound coming through your headset. To mute the television sound, press Mute on your remote control. Or, if your TV has a speaker switch on the back panel, set this switch to "OFF".

If pressing "Mute" on your remote also turns off the sound from your headset, your TV may have speaker options that you can access with your television's remote control. Please refer to your television instruction manual to locate this control and set it to "OFF".

2. **The sound from your television will only be heard through your headset and will not be heard from the television speakers.** This will happen if you connect the transmitter using the headphone adapter. If you want the sound to be heard from the headset and the television speakers, connect the transmitter to your television or cable box using the RCA stereo adapter. (refer to Pg. 6)

NOTE: If you hear static or a radio signal from your headset, it may be due to the electric current to the transmitter. Unplug the AC adapter from the electrical outlet and plug it into a surge protector. Then plug the surge protector into the electrical outlet.

If you are using your TV Listener K4 with an LCD or plasma television and are experiencing static, connect the transmitter to your cable box instead of your television. Move the transmitter as far from your television as possible while making sure it is still in the line of sight of the headset.

There is no limit to the number of headsets that can be used with a single transmitter. Please call us at **1-800-430-0222** to find out where you can purchase additional headsets.

Tips for Using Your TV Listener K4

- When connecting to your television or other audio equipment, be sure to use the audio **OUTPUT** jacks. The jacks are labeled "AUDIO OUTPUT" (not just "AUDIO"). If you are connecting to your television, these jacks are located on the back of the television. Do **NOT** use the audio jacks on the front as these are audio input jacks.
- If you are using your TV Listener K4 with an LCD or plasma television and are experiencing static, connect the transmitter to your cable box instead of your television. Follow the instructions for Stereo Installation on Pg. 6. Also, move the transmitter so that it is not directly in front of your television. Move the transmitter as far from the television as possible while making sure it is still in the line of sight of the headset.
- If you hear static or a radio signal from your headset, it may be due to the electric current to the transmitter. Unplug the AC adapter from the electrical outlet and plug it into a surge protector. Then plug the surge protector into the electrical outlet. Also, try connecting the transmitter to your cable box instead of your television.
- The TV Listener K4 does not work with optical or digital audio outputs. If your television does not have analog RCA audio outputs, use the RCA audio outputs on your cable box. If there are also no analog RCA audio outputs on your cable box, use the headphone adapter. Or you may purchase an audio converter at an electronic store.
- When charging your headset, make sure the transmitter is plugged into an electrical outlet and the headset is turned off.
- Only use rechargeable NiMH batteries. Do **NOT** use alkaline batteries or you will severely damage the headset.
- If your headset is not working, check the batteries first. Make sure the batteries are charged and properly installed.

Care & Cleaning

Disconnect the transmitter and remove all batteries from the headset before cleaning. Clean the exterior of the units with a soft, dry or slightly damp cloth. Do **NOT** spray with cleaning products or solvents. Do **NOT** submerge in water or any other liquid.

Troubleshooting

If you believe your TV Listener K4 is not operating properly, please refer to the Troubleshooting Guide below. Improper installation of the transmitter is the leading cause of malfunctions. If you are unable to solve the problem after checking the Troubleshooting Guide, please contact our Customer Service Department at **1-877-736-6760** Mon-Fri from 9:00AM-5:00PM Eastern Time, or e-mail us at customerservice@unisar.com.

PROBLEM	POSSIBLE CAUSE	SOLUTION
Four red lights on transmitter do not illuminate.	AC adapter is not connected or electrical outlet is out of order.	Check connection or change to another electrical outlet.
	There is no audio output from audio device.	Check sound on your TV or other audio device.
Power indicator light on headset does not illuminate.	Batteries in headset are low or improperly installed.	Check or charge batteries. Use only AAA rechargeable batteries. Do NOT use alkaline batteries.
	Headset is not turned on.	Make sure the power switch is in the "ON" position.
No sound from headset.	Transmitter and/or headset are not on.	Make sure both units are on.
	Batteries in headset are low or improperly installed.	Check or charge batteries. Use only AAA rechargeable batteries. Do NOT use alkaline batteries.
	Headset is not in line of sight of transmitter.	Make sure nothing is obstructing the transmission of the infrared light rays from transmitter to headset.
	No audio output from audio device.	Check sound on your audio device.
	Transmitter is not properly plugged into audio device.	Check connections. Make sure adapter is plugged into the audio output jack (not input).
	Headset and transmitter are too far apart.	Bring them closer together. Maximum range is 33 feet (10 meters).

Sound from headset is low or poor.

Volume on headset is too low.	Increase volume level.
Headset and transmitter are too far apart.	Bring them closer together. Maximum range is 33 feet (10 meters).
Headset is not in line of sight of transmitter.	Make sure nothing is obstructing the transmission of the infrared light rays from transmitter to headset.
Transmitter is not placed properly.	Place transmitter at or close to your eye level when you are sitting.
Infrared sensors in headset are obstructed.	Make sure infrared sensors are not being blocked by clothing, hair, hands, etc.
Interference from sunlight.	Position unit away from direct sunlight.
Batteries in headset are low.	Check or charge batteries. Use only AAA rechargeable batteries. Do NOT use alkaline batteries.
Electric current to transmitter is causing static.	Plug AC adapter into surge protector and then plug surge protector into electrical outlet.
LCD or plasma television is causing static.	Plug the transmitter into your cable box instead of your television. Move transmitter so it is not directly in front of the TV.
Headset is not in line of sight of transmitter.	Make sure there isn't anything obstructing the transmission of the infrared light rays from transmitter to headset.
Electric current to transmitter is causing the radio signal.	Plug AC adapter into surge protector and then plug surge protector into electrical outlet.
Connection to your television is causing the radio signal.	Plug the transmitter into your cable box instead of your television.

There is static coming from headset.

There is a radio signal coming from headset.

There is sound coming from only one side of the headset.

Audio cord is not properly plugged into audio device.

Check connections. Make sure you have not connected the stereo adapter into a mono output.

Volume on one side of headset is too low.

Increase volume level.

Technical Specifications

Power Supply > 12V, 200mA AC adapter for transmitter
2 AAA NiMH rechargeable batteries for headset

Maximum Operating Distance > 33 feet (10 meters)

Battery Life > Approximately 10 hours

Transmission Mode > Infrared 850 nanometers

Frequency Modulation > 2.3 MHz

Frequency Response > 40 Hz–12000 Hz

Distortion > < 2% THD

Signal to Noise Ratio > 40 dB

Maximum Amplification > 125-130 dB

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

This equipment has been tested and found to comply with the limits for Class B Digital Device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Any changes or modifications not expressly approved by the party responsible for compliance may void the user's authority to operate the equipment.

Limited Product Warranty

Unisar Inc., the distributor of the TV Listener® K4 (Model #TV865), warrants to the original owner of this product that this product is free from defects in materials and workmanship for a period of six (6) months from the date of the original purchase. If this product fails to function properly when used as directed under normal conditions within the six-month period, Unisar Inc. will either repair or replace the product, at its option, free of charge, when it is returned in its original packaging, postage prepaid, with proof of purchase. Unisar Inc. will bear the cost of repairing or replacing the product and shipping it back to you within the United States. Repair or replacement of the defective product is subject to verification of the malfunction or defect when delivered to:

Unisar® Service Center
51-02 21 Street, Floor 7
Long Island City, NY 11101

NOTE: The warrantor expressly disclaims liability for incidental, consequential, or special damages of any nature. This warranty does not cover: 1) A defect that has been caused by repair or modification by any person other than a Unisar technician; 2) A defect caused by negligence, accident, or misuse of the product; 3) A defect caused by incorrect installation of and/or defective batteries; 4) Battery life; 5) A defect caused by the use of any adapter not provided by the warrantor; or 6) Damage resulting from improper operation.

LEGAL RIGHTS VARY FROM STATE TO STATE.

NOTE: This warranty is valid only for products sold in the United States. If you have purchased this product outside of the United States, please see your local retailer or distributor for service.

If you have any questions, please contact our Customer Service Department at **1-877-736-6760** Mon-Fri from 9:00AM-5:00PM Eastern Time, or e-mail us at **customerservice@unisar.com**.

Also Available

TV Listener K4 Extra Headset

Model #TV865-HS

Allows multiple people to use your TV Listener K4 at the same time. No limit to the number of headsets that can be used with a single transmitter. Volume control on the headset allows each user to adjust the volume to his/her own needs.

Portable miniFRIDGE

Model #MF104

Keeps your food and beverages cold or hot. Compact and lightweight. Convenient for nursery, bedroom, office, or car. Holds baby bottles, snacks, or even six 12-ounce cans. Removable shelf. Includes home AC adapter and car adapter.

AlwaysClean™ Pacifier

Model #BR120

Closes when dropped!™ Unique shield automatically closes when the paciifer is dropped keeping it sanitary for your baby. Silicone orthodontic nipple. 0-6 months.

kickTrak™

Prenatal Kick & Movement Counter

Model #BK100

Count your way to a healthy baby!™ Your unborn baby's movement is an excellent indicator of her well-being. Makes counting and keeping a record of these movements easy and convenient with just a simple click of a button.

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